



FICHA TÉCNICA Y
DOCUMENTACIÓN

MASCARILLA QUIRÚRGICA IIR

FÁBRICA
MY TICARET VE MEDIKAL A.S

FARMA *inca*

NOMBRE DEL PRODUCTO:**Mascarilla quirúrgica tipo IIR.**

BASIC UDI DI- 868227994MSKXT

Tamaño adulto e infantil.

Desechable.

De un solo uso.

Exentas de látex.

No Estéril.

DISTRIBUIDOR:**INCA PRODUCT S.L.**

B-05115159

Calle Reyes Católicos, 8 2ºB

05001 Ávila

info@incaproduct.com

Telf.: 920269774

FABRICANTE:**MY TICARET VE MEDIKAL A.S**

Ömerli Mahallesi General Sükrü Koraltı

Caddesi n° 33 Arnavutköy

Istanbul (Turkey)

CERTIFICACIONES y DOCUMENTACIÓN ADJUNTA:

-Declaración de conformidad UE conforme Reglamento de producto sanitario 2017/745 y estándar EN 14683:2019+AC:2019

-Certificado de gestión de calidad de productos sanitario, ISO 13485:2016

- Certificado de gestión de calidad, ISO 9001:2015+AC

-Test report emitido por EKOTEKS.

-Ficha técnica My Medikal, con la siguiente documentación:

Información general de la mascarilla

Parametros de estándar EN 14683:2019+AC:2019

Test report de irritabilidad ISO 10993-10

Test report de sensibilidad ISO 10993-10

Test / inspection report por EROLAB LABORATORY SERVICES.

Ekoteks



TECHNICAL SHEET

Prepared for:

INCA PRODUCTS S.L.
C/Reyes Católicos, 8. 2º B
05001 Ávila España
Tel (+34) 616 79 40 13

Manufacturer:

My Ticaret ve Medikal A.Ş.
Ömerli Mah. General Şükrü Koraltı Cad.
No:33 Arnavutköy, İstanbul

Brand: INCA

Description: Non-sterile, Surgical Mask Type IIR (Adult and Kids)

DEFINITION : INCA 3-ply surgical mask (facial mask for medical use), can be fit according to each face measures and shapes, flexible and can be used without disturbing the soft structure. Air permeable and lets breathing easily. Non-irritating. Provides protection against bacteria.

QUALITY SYSTEM OF MY TICARET VE MEDIKAL :

MY TICARET VE MEDIKAL has been manufacturing its products with the quality systems as given below;


- EN ISO 13485:2016 Quality Assurance System-Medical Devices
- EN ISO 9001:2015+AC Quality Assurance System

PRODUCT IDENTIFICATION AND RECOMMENDED USE : The product is made by non-woven fabric. The product's composition is polypropylene and does not includes latex. The product is breathable and has no special personal on environmental hazards. The product is made automatically in hygienic conditions. The product prevents the potential reactions between all kind of liquids and particles, microorganisms.

SHELF LIFE : 5 Years

DIMENSION

- Adult Surgical Mask

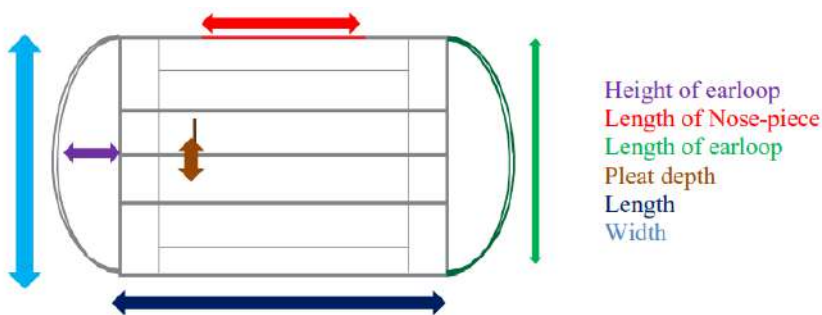

| | | | |
|------------------|-----------------------------|---------------|--------|
| Dimension | Body Size | Length | 175 mm |
| | | Width | 95 mm |
| | Height of Earloop | | 70 mm |
| | Lenght of Earloop | | 160 mm |
| | Pleat Depth | | 10 mm |
| | Length of Nose-piece | | 90 mm |

- Kids Masks

| | | | |
|-------------------------------|-----------------------------|-------------------------|--------|
| Dimension (Infant) | Body Size | Length (+/-1 cm) | 145 mm |
| | | Width (+/-1 cm) | 95 mm |
| | Height of Earloop | | 70 mm |
| | Length of Earloop | | 160 mm |
| | Pleat Depth | | 10 mm |
| | Length of Nose-piece | | 60 mm |

| Characteristic | Specification | | |
|-----------------------|----------------------|-----------------------|--------------------|
| Materials | Mask Body | Outer Material | Spunbond 25 -30 gr |
| | | Filter Layer | Meltblown 25 gr |
| | | Inner Material | Spunbond 20-25 gr |

Reference/Catalog No. and Barcodes

| Description | Quantity /box | Catalog # | Inner Box EAN Barcode | Inner box UDI Barcode |
|---|----------------------|------------------|------------------------------|--|
| Adult Masks, Surgical Type IIR (Black) | 50pcs | MM.NS.LM.510 | 8683020020603 | (01)08683020020603(11)MMYY(10)XXXXXXXXXX |
| Adult Masks, Surgical Type IIR (Blue) | 10pcs | MM.NS.LM.511 | 8683020020610 | (01)08683020020610(11)MMYY(10)XXXXXXXXXX |
| Kids Mask, Surgical Type IIR (Blue) | 50pcs | MM.NS.LM.512 | 8683020020627 | (01)08683020020627(11)MMYY(10)XXXXXXXXXX |
| Kids Mask, Surgical Type IIR (Blue) | 10pcs | MM.NS.LM.513 | 8683020020634 | (01)08683020020634(11)MMYY(10)XXXXXXXXXX |
| Adult Masks, Surgical Type IIR (Blue) | 50pcs | MM.NS.LM.514 | 8683020020641 | (01)08683020020641(11)MMYY(10)XXXXXXXXXX |
| Kids Mask, Surgical Type IIR (Pink) | 10pcs | MM.NS.LM.558 | 8683020020665 | (01)08683020020665(11)MMYY(10)XXXXXXXXXX |
| Kids Mask, Surgical Type IIR (Black) | 10pcs | MM.NS.LM.559 | 8683020020672 | (01)08683020020672(11)MMYY(10)XXXXXXXXXX |
| Kids Mask, Surgical Type IIR (Octopus Printed) | 10pcs | MM.NS.LM.560 | 8683020020689 | (01)08683020020689(11)MMYY(10)XXXXXXXXXX |
| Kids Mask, Surgical Type IIR (Panda Printed) | 10pcs | MM.NS.LM.561 | 8683020020696 | (01)08683020020696(11)MMYY(10)XXXXXXXXXX |

MY TİCARET VE MEDİKAL A.Ş.

Ömerli Mah. General Şükrü Koraltı Cad. No:33 Arnavutköy / İstanbul

Tel: 0212 438 20 64

Fax: 0212 438 20 65

info@mymedikal.com.tr

www.mymedikal.com.tr



MY Medikal

Loading Information (Mega Truck)

| Brand | Description | Pcs per box | No. of box per carton | No. of carton per pallet | Max Qty (ctn) | Quantity (boxes) |
|-------|--------------------------------|-------------|-----------------------|--------------------------|---------------|------------------|
| INCA | Adult Masks, Surgical Type IIR | 50pcs | 60 | 16 | 660 | 39,600 |
| INCA | Adult Masks, Surgical Type IIR | 10pcs | 210 | 16 | 660 | 138,600 |
| INCA | Kids Mask, Surgical Type IIR | 50pcs | 60 | 20 | 770 | 46,200 |
| INCA | Kids Mask, Surgical Type IIR | 10pcs | 180 | 20 | 770 | 138,600 |

| Description | Pcs per box | No. of box per carton | Inner box Dimension | Outer carton Dimension |
|--------------------------------|-------------|-----------------------|--------------------------|------------------------|
| Adult Masks, Surgical Type IIR | 50pcs | 60 | 18,5cm * 10cm * 10,5cm | 60cm * 40cm * 56cm |
| Adult Masks, Surgical Type IIR | 10pcs | 210 | 10.5cm * 20.0cm * 2.68cm | 60cm * 40cm * 56cm |
| Kids Mask, Surgical Type IIR | 50pcs | 60 | 15,4cm * 10cm * 10,4cm | 60cm * 33cm * 56cm |
| Kids Mask, Surgical Type IIR | 10pcs | 180 | 10,5cm * 17cm * 2,925cm | 60cm * 33cm * 56cm |

*Remarks: Maximum quantity refers to unpalletized goods and loaded to a mega truck.
Mega Truck dimension: 240cm – width x 280cm – height x 1360cm – length*

For loading with pallets, information:

| | |
|---|-------------|
| Pallet Dimension for adult masks (cm): | 80*120*240 |
| Pallet Dimension for kids masks (cm): | 100*120*240 |
| Pallet Weight (without masks): | 13 kgs |
| No. of Pallets per Truck(for adult masks): | 33 |
| No. of Pallets per Truck(for kids masks): | 26 |

MY TİCARET VE MEDİKAL A.Ş.

Ömerli Mah. General Şükrü Koraltı Cad. No:33 Arnavutköy / İstanbul

Tel: 0212 438 20 64

Fax: 0212 438 20 65

info@mymedikal.com.tr

www.mymedikal.com.tr



MANAGEMENT SYSTEM
ISO / IEC 17021-1:2015
NAC-002-MS

Certificate of Registration

This is to certify that

Quality Management System

of

MY TİCARET VE MEDİKAL ANONİM ŞİRKETİ

**ÖMERLİ MAHALLESİ GENERAL ŞÜKRÜ KORALTI CAD. NO: 33
ARNAVUTKÖY - İSTANBUL / TÜRKİYE**

complies with requirements of

ISO 9001:2015

This certificate is valid concerning all activities related to;

PRODUCTION AND SALE OF POWDERED/POWDER-FREE LATEX EXAMINATION GLOVES, POWDER-FREE NITRILE EXAMINATION AND PROTECTIVE GLOVES, POWDERED/POWDER-FREE STERILE SURGICAL GLOVES, POWDERED/POWDER-FREE VINYL EXAMINATION GLOVES, SURGICAL MASK, FILTERING HALF/PROTECTIVE FACE MASK, STERILE & NON-STERILE GAUZE SWABS/ LAP SPONGES/ COTTON PADS/ GAUZE, STERILE SYRINGE/NEEDLE

PUDRALI/PUDRASIZ LATEKS MUAYENE ELDİVENİ, NİTRİL PUDRASIZ MUAYENE VE KORUYUCU ELDİVEN, PUDRALI VE PUDRASIZ STERİL CERRAHİ ELDİVEN, VİNİL PUDRALI/PUDRASIZ MUAYENE ELDİVENİ, CERRAHİ MASKE, FİLTRELİ KORUYUCU YÜZ MASKESİ, STERİL & NON-STERİL GAZ KOMPRES/BATIN KOMPRES/PAMUKLU PED/GAZLI BEZ, STERİL ŞİRINGA/İĞNE ÜRETİMİ VE SATIŞI

ISO 01 940 1179
Certificate No.

Jun. 2, 2021
Date of this Certificate

Jun. 4, 2022
Certification Expiry Date

Apr. 5, 2021
Date of Audit

Jun. 5, 2020
Date of Registration


Managing Director / Director



Medicert Uluslararası Ürün Ve Sistem Belgelendirme Ltd. Şti.
Tersane Mah. Cemal Gürsel Cad. No:11/3 Halide Hnm. Apt. Karşıyaka / İzmir
Tel: 0232 327 33 44 www.medicert.com.tr info@medicert.com.tr

This certificate is only valid if it is available/valid on Medicert website at www.medicert.com.tr

This certificate of Registration remains the property of Medicert Certificate Ltd and shall be returned immediately upon request
* In Case if Surveillance Audit is not allowed to be conducted on or before the specified date; the Certificate shall be Suspended/Withdrawn.





Certificate of Registration

This is to certify that

**Quality Management System
for Medical Devices**

of

MY TİCARET VE MEDİKAL ANONİM ŞİRKETİ

**ÖMERLİ MAHALLESİ GENERAL ŞÜKRÜ KORALTI CAD. NO: 33
ARNAVUTKÖY - İSTANBUL / TÜRKİYE**

complies with requirements of

ISO 13485:2016

This certificate is valid concerning all activities related to;

PRODUCTION AND SALE OF POWDERED/POWDER-FREE LATEX EXAMINATION GLOVES, POWDER-FREE NITRILE EXAMINATION AND PROTECTIVE GLOVES, POWDERED/POWDER-FREE STERILE SURGICAL GLOVES, POWDERED/POWDER-FREE VINYL EXAMINATION GLOVES, SURGICAL MASKS, PROTECTIVE FACE MASK WITH FILTER

PUDRALI/PUDRASIZ LATEKS MUAYENE ELDİVENİ, NİTRİL PUDRASIZ MUAYENE VE KORUYUCU ELDİVEN, PUDRALI VE PUDRASIZ STERİL CERRAHİ ELDİVEN, VİNİL PUDRALI/PUDRASIZ MUAYENE ELDİVENİ, CERRAHİ MASKE, FİLTRELİ KORUYUCU YÜZ MASKESİ ÜRETİMİ VE SATIŞI

ISO 02 836 1179
Certificate No.

Mar. 8, 2021
Date of this Certificate

Feb. 25, 2022
Certification Expiry Date

Feb. 8, 2021
Date of Audit

Feb. 26, 2020
Date of Registration


Managing Director / Director



Medicert Uluslararası Ürün Ve Sistem Belgelendirme Ltd. Şti.
Tersane Mah. Cemal Gürsel Cad. No:11/3 Halide Hnm. Apt. Karşıyaka / İzmir
Tel: 0232 327 33 44 www.medicert.com.tr info@medicert.com.tr

This certificate is only valid if it is available/valid on Medicert website at www.medicert.com.tr

This certificate of Registration remains the property of Medicert Certificate Ltd and shall be returned immediately upon request
* In Case if Surveillance Audit is not allowed to be conducted on or before the specified date; the Certificate shall be Suspended/Withdrawn.



**DECLARATION OF CONFORMITY TO
MEDICAL DEVICE REGULATION 2017/745**

Manufacturer : My Ticaret ve Medikal A.Ş
Manufacturer's Address : Ömerli Mahallesi General Şükrü Koraltı
Caddesi No: 33 Arnavutköy/İstanbul/Turkey

This declaration of conformity covers the Class 1 medical device products mentioned below.

Device Description:

| Brand | Device Name | Catalog /Ref No |
|--------------|--|------------------------|
| INCA | Non Sterile Adult Surgical Masks, 3 Ply Type IIR (Black) 50 pcs | MM.NS.LM.510 |
| INCA | Non Sterile Adult Surgical Masks, 3 Ply Type IIR (Blue) 10 pcs | MM.NS.LM.511 |
| INCA | Non-Sterile Kids Surgical Mask, 3 Ply Type IIR (Blue) 50 pcs | MM.NS.LM.512 |
| INCA | Non- Sterile Kids Surgical Mask, 3 Ply Type IIR (Blue) 10 pcs | MM.NS.LM.513 |
| INCA | Non-sterile Adult Surgical Masks, 3 Ply Type IIR (Blue) 50 pcs | MM.NS.LM.514 |
| INCA | Non-sterile Kids Surgical Mask, 3 Ply Type IIR (Pink) 10 pcs | MM.NS.LM.558 |
| INCA | Non-sterile Kids Surgical Mask, 3 Ply Type IIR (Black) 10 pcs | MM.NS.LM.559 |
| INCA | Non-sterile Kids Surgical Mask, 3 Ply Type IIR (Octopus Printed) 10 pcs | MM.NS.LM.560 |
| INCA | Non-sterile Kids Surgical Mask, 3 Ply Type IIR (Panda Printed) 10 pcs | MM.NS.LM.561 |

A surgical mask, also known as a medical face mask, is a personal protective equipment worn by healthcare professionals during medical procedures. It is a medical device that, when properly fitted, prevents the transmission of airborne infections between patients and/or treatment personnel by preventing pathogens (primarily bacteria and viruses) emitted in respiratory droplets and aerosols from entering the user's mouth and nose.

Device Classification : Medical Device Class I, according to MD Regulation (EU)2017/745.
Conformity Assesment Procedure : Annex II&III Medical Device Regulation (EU)2017/745
MDR Conformity Route: : Self Declaration
Basic UDI-DI : 868302002MSKMH

We, My Ticaret ve Medikal A.Ş. herewith declared that above mentioned device:

- Is in conformity with The Regulation (EU) 2017/745 of The European Parliament and of The Council of medical devices.

This EU declaration of conformity is issued under the sole responsibility of the manufacturer, My Ticaret ve Medikal A.Ş.

Applied standards: EN 14683:2019+AC:2019

Verification Certificates Authorized Signatory

: Quality Management System EN
13485:2016

EU DoC Issuance date

: 30.11.2021

Signed For and Behalf of My Ticaret ve Medikal A.Ş.

Name: MURAT YILDIZ

Designation: General Manager

MY TICARET VE
MEDİKAL A.Ş. ŞİRKETİ
Ömerli Mah. General Şükrü Koraltı Cad
No 33 Arnavutköy/İSTANBUL
Buyukörsençe V.D. 626 040 4605
Tel: 0212 438 20 64 Fax: 0212 438 20 65
www.mymedikal.com.tr



**EKOTEKS LABORATUVAR ve GÖZETİM
HİZMETLERİ A. Ş.**
Esenyurt Firuzköy Bulvarı No:29 34325 Avcılar
İstanbul/ TÜRKİYE



TEST REPORT
DENEY RAPORU

| |
|------------------|
| AB-0583-T |
| 20031582- ing |
| 09-20 |

Customer name: MY TİCARET VE MEDİKAL A.Ş.
Address: Ömerli Mahallesi General Şükrü Koraltı Caddesi No:33 ARNAVUTKÖY/
İSTANBUL
Buyer name: -
Contact Person: Z.MELEK ÖZ BOLAT
Order No: SAMPLE 1
Article No: MM.NS.LM.01 MUMU SURGICAL MASK
Name and identity of test item: Blue non-woven mask.
The date of receipt of test item: 01.09.2020
Re-submitted/re-confirmation date: -
Date of test: 01.09.2020-10.09.2020
Remarks: -
Sampling: The results given in this report belong to the received sample by vendor.
End-Use: -
Care Label: Not specified.
Number of pages of the report: 5

The Turkish Accreditation Agency (TURKAK) is signatory to the multilateral agreements of the European co-operation for the Accreditation (EA) and of the International Laboratory Accreditation (ILAC) for the Mutual recognition of test reports.

EKOTEKS LABORATUVAR ve GÖZETİM HİZMETLERİ A.Ş. accredited by TÜRKAK under registration number [AB-0583-T] for ISO 17025:2017 as test laboratory.

The test and/or measurement results, the uncertainties (if applicable) with confidence probability and test methods are given on the following pages which are part of this report.

Seal

Date
11.09.2020

Customer Representative
Özlem ULUS

Head of Testing Laboratory
Sevim A. RAZAN
11.09.2020

This report shall not be reproduced other than in full except with the permission of the laboratory.
Testing reports without signature and seal are not valid.

EKOTEKS LABORATUVAR ve GÖZETİM
HİZMETLERİ A.Ş.


AB-0583-T

20031582-
ing

09-20

| REQUIRED TESTS | RESULT | COMMENTS |
|--|--------|----------|
| MICROBIOLOGICAL TESTS | | |
| Bacterial Filtration Efficiency-BFE | P | TYPE IIR |
| Microbial Cleanliness(Bioburden) | P | |
| Splash Resistance | P | |
| PHYSICAL PROPERTIES | | |
| Breathability(Differential Pressure) | P | |
| P: Pass F: Fail R: Refer to retailer technologist. Tests results were evaluated according to EN 14683:2019+AC :2019 Tablo 1 limit values. | | |

REMARK: Original samples are kept for 3 months and all technical records are kept for 5 years unless otherwise specified. If requested, measurement uncertainty will be reported. But unless otherwise specified, measurement uncertainty is not considered while stating compliance with specification or limit values. The reported uncertainty is based on a standard uncertainty multiplied by a coverage factor $k=2$, providing a level of confidence of approximately 95 %. Tests marked (*) in this report are not included in the accreditation schedule.

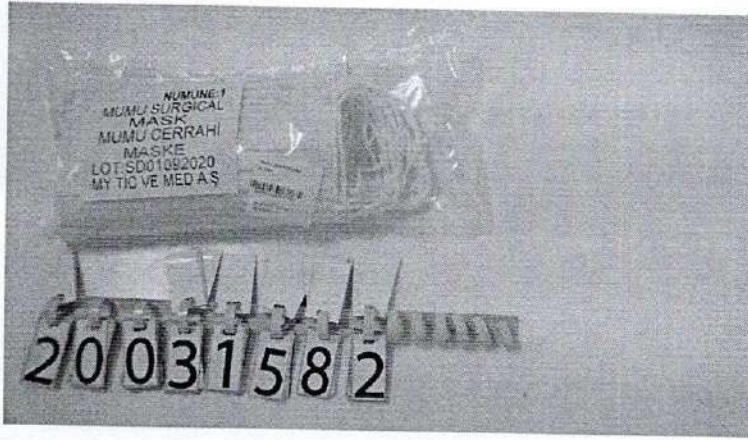

Gen.f136-2/03

This report shall not be reproduced other than in full except with the permission of the laboratory.
Testing reports without signature and seal are not valid.

AB-0583-T

20031582-
ing

09-20

TEST RESULTS

BACTERIAL FILTRATION EFFICIENCY (BFE)

Test Metod: (Bacterial Filtration Efficiency Testing –BFE /Ref: EN 14683:2019+AC:2019 Medical Face Masks, Requirements and Test Methods

A specimen of the mask material is clamped between an impactor and an aerosol chamber. An aerosol of *Staphylococcus aureus* is introduced into the aerosol chamber and drawn through the mask material and the impactor under vacuum. The bacterial filtration efficiency of the mask is given by the number of colony forming units passing through the medical face mask material expressed as a percentage of the number of colony forming units present in the challenge aerosol.

| | |
|--|--|
| Test Flow Rate | 28,3 L/min |
| Test Flow Time | 2 minute |
| Sample Sizes | 5 pieces mask |
| Microorganism | <i>Staphylococcus aureus</i> ATCC 6538 |
| Bacterial concentration (cfu/ ml) | 5×10^5 cfu/ ml |
| incubation conditions | 24 hour, 35°C ± 2°C |
| Positive control sample average of number of Bacteria (C) | $1,93 \times 10^3$ cfu/ ml |

| RESULTS | | | |
|-----------------------|---|--|-------------------------------|
| Number of Test Sample | Test Sample (T) Number of Bacteria (cfu/ml) | Bacterial Filtration Efficiency (% B) | Requirement BFE (%) |
| 1 | 38 | 98,0% | Type I ≥95 Type II ≥98 |
| 2 | 33 | 98,3% | |
| 3 | 29 | 98,5% | |
| 4 | 24 | 98,8% | |
| 5 | 31 | 98,4% | |

cfu: Colony-forming unit

$$B = (C - T) / C \times 100$$

%B: Bacterial Filtration Efficiency

C: is the mean of the total plate counts for the two positive control runs

T: is the total plate count for the test specimen

AB-0583-T

20031582-
ing

09-20

TEST RESULTS

MICROBIAL CLEANLINESS (Bioburden)

Test Metod: EN ISO 11737-1:2018

The sample is put in extraciton liquid after shaking well, inoculated on the agar.
After incubation at 30 ± 1 ° C for 72 hours, growth microorganisms are counted on the agar.

| | <u>RESULTS</u> | <u>REQUIREMENT</u> |
|--------------------------------------|----------------|--------------------|
| Microbial cleanliness (cfu/g) | 16 cfu/g | ≤ 30 cfu/g |

*cfu= Colony forming unit.

SPLASH RESİSTANCE

Test Metod: EN 14683:2019+AC :2019 (Clause 5.2.4) the resistance of the medical face mask to penetration of splashes of liquid shall conform to the minimum value given for Type IIR in Table 1

ISO 22609 :2004 Clothing for protection against infectious agents — Medical face masks — Test method for resistance against penetration by synthetic blood (fixed volume, horizontally projected)

Test Condition (21 ± 5) °C ve (85 ± 5) % relative humidity, 4 hrs

6 different samples were taken

| | <u>SPLASH RESİSTANCE PRESSURE (kPa)</u> | <u>RESULTS</u> | <u>REQUIREMENT</u> |
|-----------------------|---|----------------|--------------------|
| 1 | >21.3 kPa | PASS | ≥ 16 kPa |
| 2 | >21.3 kPa | PASS | |
| 3 | >21.3 kPa | PASS | |
| 4 | >21.3 kPa | PASS | |
| 5 | >21.3 kPa | PASS | |
| 6 | >21.3 kPa | PASS | |
| Average Result | >21.3 kPa | PASS | |

EKOTEKS LABORATUVAR ve GÖZETİM
HİZMETLERİ A.Ş.

| |
|------------------|
| AB-0583-T |
| 20031582- ing |
| 09-20 |

TEST SONUÇLARI

BREATHABILITY (Differential Pressure)

Test Method: EN 14683:2019+AC :2019 (TS EN 14683+AC:2019) Annex-C

Test Condition (21 ± 5) °C ve (85 ± 5) % relative humidity, 4 hrs

Test area is 25 mm in diameter , 5 different sample was taken

Adjusted airflow is 8 l/min. The differential pressure is read directly using a differential pressure manometer .

| SAMPLE | DIFFERENTIAL PRESSURE RESULT | REQUIREMENT |
|-----------------------|---------------------------------|-------------------------|
| 1 | 24.3 Pa/cm ² | < 60 Pa/cm ² |
| 2 | 23,1 Pa/cm ² | |
| 3 | 24,8 Pa/cm ² | |
| 4 | 29,0 Pa/cm ² | |
| 5 | 22,7 Pa/cm ² | |
| | | |
| Average Result | 24,7 Pa/cm ² | |

TEST REPORT
DENEY RAPORU

AB-0583-T

21034818

12-21

Customer name: MY TİCARET VE MEDİKAL A.Ş.
Address: Ömerli Mahallesi General Şükrü Koraltı Caddesi No:33 ARNAVUTKÖY/
İSTANBUL
Buyer name: -
Contact Person: HİLAL YİĞİT
Order No: -
Article No: Surgical Kids Mask / MM.NS.LM.29
Name and identity of test item: Printed white non-woven medical mask
The date of receipt of test item: 17.11.2021
Re-submitted/re-confirmation date: 24.11.2021
Date of test: 24.11.2021-02.12.2021
Remarks: -
Sampling: The results given in this report belong to the received sample by vendor.
End-Use: -
Care Label: Not specified.
Number of pages of the report: 3

The Turkish Accreditation Agency (TURKAK) is signatory to the multilateral agreements of the European co-operation for the Accreditation (EA) and of the International Laboratory Accreditation (ILAC) for the Mutual recognition of test reports.

EKOTEKS LABORATUVAR ve GÖZETİM HİZMETLERİ A.Ş. accredited by TÜRKAK under registration number [AB-0583-T] for ISO 17025:2017 as test laboratory.

The test and/or measurement results, the uncertainties (if applicable) with confidence probability and test methods are given on the following pages which are part of this report.



Date
02.12.2021

Customer Representative
Özlem ULUS

Head of Testing Laboratory
Sevim A. RAZAK
02.12.2021

This report shall not be reproduced other than in full except with the permission of the laboratory.
Testing reports without signature and seal are not valid.

EKOTEKS LABORATUVAR ve GÖZETİM
HİZMETLERİ A.Ş.


AB-0583-T

21034818

12-21

| REQUIRED TESTS | EVALUATION | COMMENTS |
|--|------------|----------|
| MICROBIOLOGICAL TESTS | | |
| Bacterial Filtration Efficiency-BFE | P | |
| P: Pass F: Fail R: Refer to retailer technologist. Tests results were evaluated according to EN 14683:2019+AC :2019 Tablo 1 limit values. | | |

REMARK: Original samples are kept for 3 months and all technical records are kept for 5 years unless otherwise specified. If requested, measurement uncertainty will be reported. But unless otherwise specified, measurement uncertainty is not considered while stating compliance with specification or limit values. The reported uncertainty is based on a standard uncertainty multiplied by a coverage factor $k=2$, providing a level of confidence of approximately 95 %. The declaration of conformity was given in accordance with the Simple Acceptance Decision Rule. Tests marked (*) in this report are not included in the accreditation schedule.

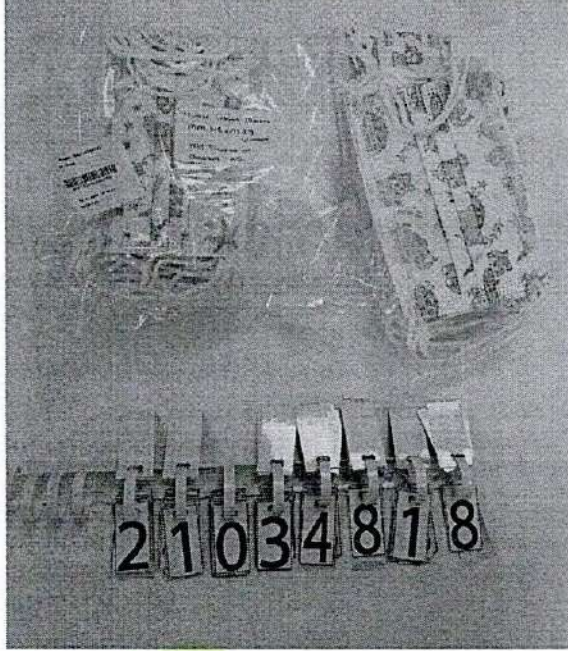

Gen.f136-2/03

This report shall not be reproduced other than in full except with the permission of the laboratory.
Testing reports without signature and seal are not valid.

**EKOTEKS LABORATUVAR ve GÖZETİM
HİZMETLERİ A.Ş.**

AB-0583-T

21034818

12-21

TEST RESULTS

BACTERIAL FILTRATION EFFICIENCY (BFE)

Test Metod: (Bacterial Filtration Efficiency Testing –BFE /Ref: EN 14683:2019+AC:2019 Medical Face Masks, Requirements and Test Methods

A specimen of the mask material is clamped between a impactor and an aerosol chamber. An aerosol of *Staphylococcus aureus* is introduced into the aerosol chamber and drawn through the mask material and the impactor under vacuum. The bacterial filtration efficiency of the mask is given by the number of colony forming units passing through the medical face mask material expressed as a percentage of the number of colony forming units present in the challenge aerosol.

| | |
|--|--|
| Test Flow Rate | 28,3 L/min |
| Test Flow Time | 2 minute |
| Sample Sizes | 5 pieces mask |
| Microorganism | <i>Staphylococcus aureus</i> ATCC 6538 |
| Bacterial concentration (cfu/ ml) | 5×10^5 cfu/ ml |
| incubation conditions | 24 hour, 35°C ± 2°C |
| Positive control sample average of number of Bacteria (C) | $3,0 \times 10^3$ cfu/ ml |

| RESULTS | | | |
|-----------------------|---|--|-------------------------------|
| Number of Test Sample | Test Sample (T) Number of Bacteria (cfu/ml) | Bacterial Filtration Efficiency (% B) | Requirement BFE (%) |
| 1 | 42 | 98,6% | Type I ≥95 Type II ≥98 |
| 2 | 44 | 98,5% | |
| 3 | 43 | 98,6% | |
| 4 | 45 | 98,5% | |
| 5 | 48 | 98,4% | |

cfu: Colony-forming unit

$$B = (C - T) / C \times 100$$

%B: Bacterial Filtration Efficiency

C: is the mean of the total plate counts for the two positive control runs

T: is the total plate count for the test specimen

Total Uncertainty: 0,40



Test Result : PASS / Level 2

Report No : 2020080611

Applicant: MY TİCARET VE MEDİKAL A.Ş.
Ömerli mah. General Şükrü Karaltı cad. No:33 Arnavutköy/ İstanbul
Z. Melek ÖZ BOLAT

Contact Person : Z. Melek ÖZ BOLAT

Contact Telephone: 0212 438 2064

Contact e-mail: Kalite@mymedikal.com.tr

Sample Accepted on: 10.07.2020

Report Date: 06.08.2020

Total number of pages: 5 (syf)

Sample ID: Surgical Mask 3-Ply With Earloop

| | TEST | METHOD | RESULTS |
|---|---|-------------|---------|
| * | Standard Specification for Performance of Materials Used in Medical Facial Masks | ASTM F 2100 | PASS |
| * | Standard Test Method for Determining the First Efficiency of Substances Used in Medical Facial Masks to Penetration Using Particles Using Latex Spheres | ASTM F 2299 | PASS |



Seal

Customer Representative
Hasan KUTLU



Laboratory Manager
Hava Sariaydin

EUROLAB[®] (TÜRCERT TEKNİK KONTROL VE BELGELENDİRME A.Ş.)

It is prohibited to change any and all versions of this document in any manner whatsoever. In case of a conflict between the electronic version (e.g. PDF file) and the original paper version provided by EUROLAB[®], the latter will prevail.

TÜRCERT Teknik Kontrol ve Belgelendirme A.Ş. disclaim liability for any direct, indirect, consequential or incidental damages that may result from the use of the information or data, or from the inability to use the information or data contained in this document.

The contents of this report may only be transmitted to third parties in its entirety and provided with the copyright notice, prohibition to change, electronic versions' validity notice and disclaimer.

Environment

The requirements and standards apply to equipment intended for use in

| | |
|---|---|
| X | Residential (domestic) environment |
| X | Commercial and light-industrial environment |
| X | Industrial environment |
| X | Medical environment |



ASTM F 2100 : Standard Specification for Performance of Materials Used in Medical Facial Masks

This specification covers tests and requirements for the materials used in making medical face masks used in the provision of health services such as surgery and patient care.

This specification provides a classification of medical face mask material performance. Medical face mask material performance is based on bacterial filtration efficiency, differential pressure, submicron particle filtration efficiency, resistance to synthetic blood penetration and flammability tests.

Necessity :

| Characteristic | Level 1 | Level 2 | Level 3 |
|---|---------|---------|---------|
| Bacterial Filtration Efficiency (%) | % 95 | % 98 | % 98 |
| Differential Pressure mm H ₂ O / cm ² | < 4,0 | < 5,0 | < 5,0 |
| Submicron Particle Filtration Efficiency 0.1 micron (%) | % 95 | % 98 | % 98 |
| Synthetic Blood Penetration Resistance (min.pressure in mmHg) | 80 | 120 | 160 |
| Combustion Class | Class 1 | Class 1 | Class 1 |



ASTM F 2299 : Standard Test Method for Determining the First Efficiency of Substances Used in Medical Facial Masks to Penetration Using Particles Using Latex Spheres

This test method creates procedures to measure the initial particle filtration efficiency of materials used in medical face masks using monodispers aerosols. This test method uses luminous particle counts from 0.1 to 5.0 μm in size and air flow test speeds of 0.5 to 25 cm / s.

The test procedure measures the filtration efficiency by comparing the number of particles in the feed stream (upstream) with the filtrate (downstream).



Test Results

| Tests | Results | | | |
|--|---------------------|---------------------|---------------------|---------------------|
| | 1 | 2 | 3 | Avg |
| Liquid Resistance, mmHg | 116 | 116 | 116 | 116 |
| BFE | % 98,85 | % 98,82 | % 98,91 | % 98,86 |
| PFE, 0.1 micron % | % 98,01 | % 98,06 | % 98,04 | %98,03 |
| Differential Pressure, mm H ₂ O / cm ² | 2,3 | 2,2 | 2,2 | 2,23 |
| Combustion Class | Class 1 Flare up | Class 1 Flare up | Class 1 Flare up | Class 1 Flare up |

Comment : Results meet Level 2 criteria



Images Under Test

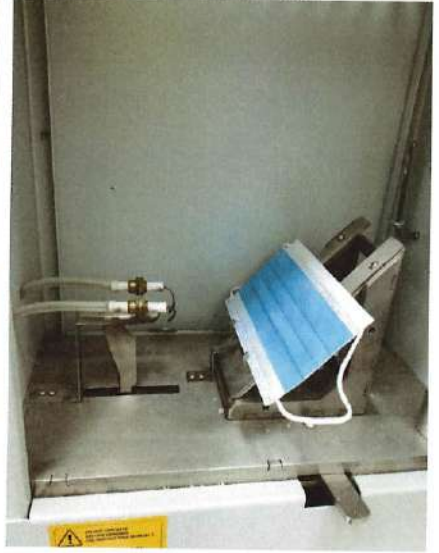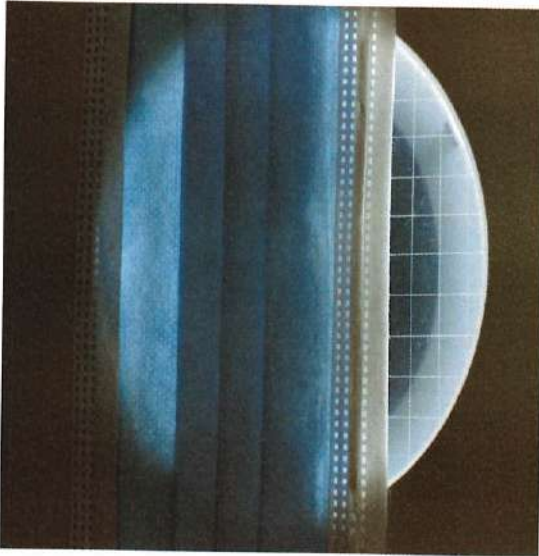

End of Report



PACKING METHOD

Ref No. MM.NS.LM.512



Ref No. MM.NS.LM.510



MY TİCARET VE MEDİKAL A.Ş.

Ömerli Mah. General Şükrü Koraltı Cad. No:33 Arnavutköy / İstanbul

Tel: 0212 438 20 64

Fax: 0212 438 20 65

info@mymedikal.com.tr

www.mymedikal.com.tr

Ref No. MM.NS.LM.513



Master carton



MY TİCARET VE MEDİKAL A.Ş.

Ömerli Mah. General Şükrü Koraltı Cad. No:33 Arnavutköy / İstanbul

Tel: 0212 438 20 64

Fax: 0212 438 20 65

info@mymedikal.com.tr

www.mymedikal.com.tr